


C-Spine WO Contrast

Contrast:N/A			
Saline Chaser: N/A			
Contrast amount: N/A			
Rate: N/A			
Delay: N/A			
Oral: None			
Display Field of View: 140mm			
Coverage:			
bottom of T1 through top of mastoids			
Images to send to PACS (# of series)			
Plane	Thickness	Kernel	Window
Axial	1mm	vendor appropriat bone kernel	osteo
Axial	3mm	vendor appropriat soft tissue kernel	mediastinum
Sagittal	2mm	vendor appropriat bone kernel	osteo
Sagittal	3mm	vendor appropriat soft tissue kernel	mediastinum
Coronal	2mm	vendor appropriat bone kernel	osteo
Comments			
If a lot of metal artifact do an additional metal reduction technique			

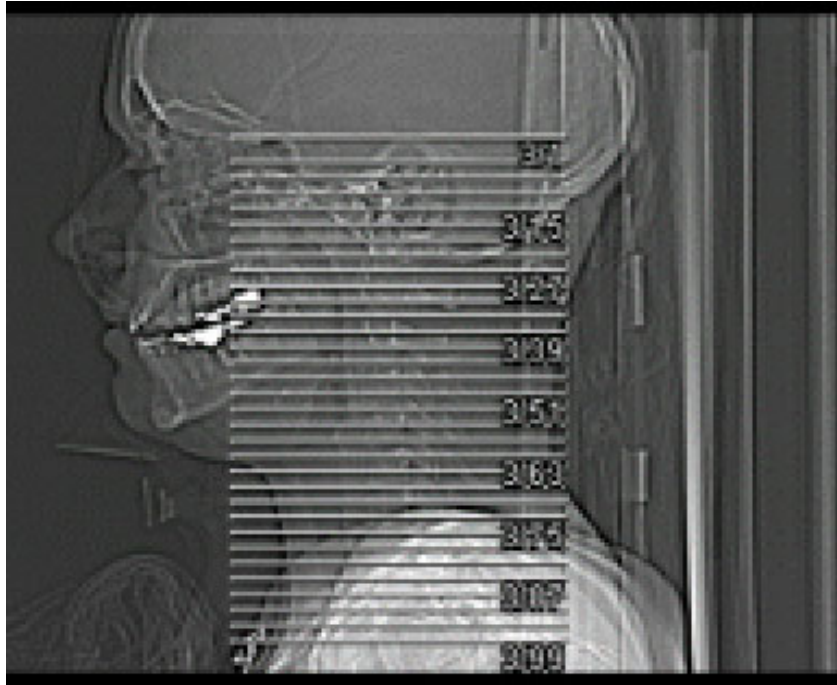

Access our Modality Protocols site for all current protocols at:

www.MIARAD.com

INTEGRITY | TEAMWORK | EXCELLENCE | SERVICE | ACCOUNTABILITY

Updated 05/2025